

LOCAL




NATIONAL



AESTHETIC




TECHNOLOGY

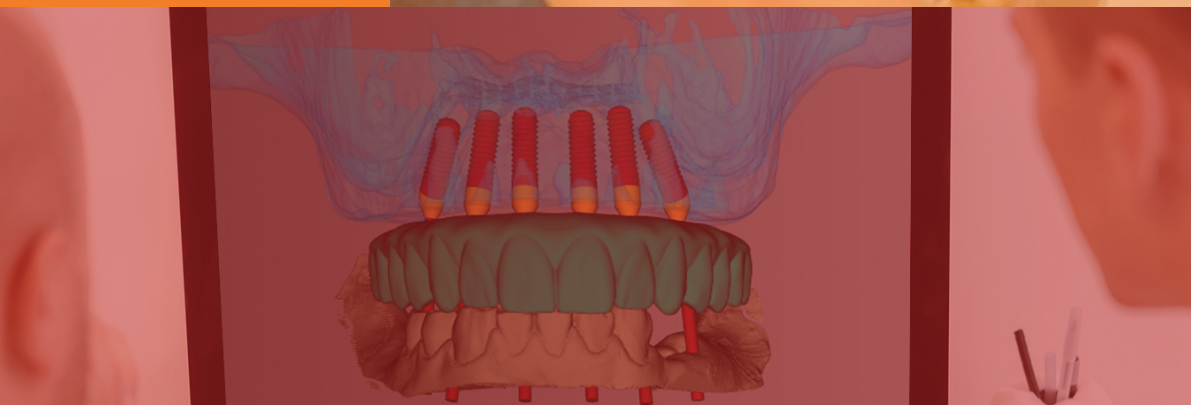

One Partner, Multiple Options

Our professionals have the unique ability to provide localized, consultative services on all of the products supported by our national network of laboratories and an unmatched array of innovative, restorative solutions, services and technologies.

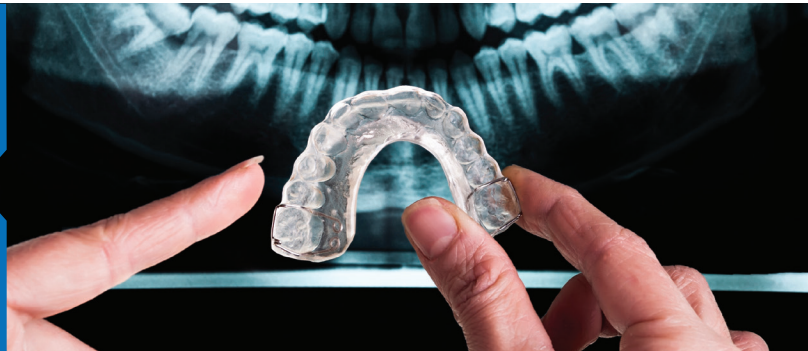

NTI^{TSS}
YOUR MIGRAINE GUARD

LAB INFORMATION

National Dentex Local Full Service Labs

NDX Affinity	2900 Blankenbaker Parkway, Suite 110	Louisville	KY	40299	502.499.7400	NDXAffinity.com
NDX Albensi	100 Colony Drive	Irwin	PA	15642	800.734.3064	NDXAlbensilab.com
NDX Bauer	1415 North Sanborn Boulevard	Mitchell	SD	57301	605.996.6521	NDXBauer.com
NDX Continental	960 W Behrend Drive, Suite 4	Phoenix	AZ	85027	623.580.7711	NDXContinental.com
NDX Dental Art	1721 N Grand River Avenue	Lansing	MI	48906	800.444.3744	NDXDentalArt.com
NDX Dexterity	310 Division Street	Farmington	MN	55024	651.463.4444	NDXDexterity.com
NDX Dodd	24 Lukens Drive	New Castle	DE	19720	302.661.6000	NDXDodd.com
NDX Florida	4520 Parkbreeze Court	Orlando	FL	32808	407.781.4725	NDXFlorida.com
NDX Green	1099 Wilburn Road	Heber Springs	AR	72543	501.362.3132	NDXGreen.com
NDX H&O	1050 Perimeter Road, Suite 101	Manchester	NH	03103	603.624.4339	NDXH-O.com
NDX Ideal	5501 Wilshire NE, Suite B	Albuquerque	NM	87113	505.265.2200	NDXIdeal.com
NDX Keller	160 Larkin Williams Industrial Court	Fenton	MO	63026	636.600.4200	NDXKeller.com
NDX Kiess Kraft	6601 S 118th Street	Omaha	NE	68137	402.391.8424	NDXKiessKraft.com
NDX Lords	1500 West Main Avenue	De Pere	WI	54115	920.347.2929	NDXLords.com
NDX Lumident	8840 Commerce Park Place, Building 4	Indianapolis	IN	46268	866.586.4336	NDXLumident.com
NDX Mallow-Tru	2650 NE Hagan Road	Lee's Summit	MO	64064	816.220.7800	NDXMallowTru.com
NDX Oral Arts Iowa	6701 Chavenelle Road	Dubuque	IA	52002	563.556.1911	NDXOralArtsIowa.com
NDX Peterman	8020 Safari Drive	Smyrna	TN	37167	615.331.1670	NDXPeterman.com
NDX Rogers	2407 Denso Drive	Athens	TN	37303	423.745.7115	NDXRogers.com
NDX Saber	803 West College Avenue	Waukesha	WI	53186	262.542.5100	NDXSaber.com
NDX Salem	3873 Rocky River Drive	Cleveland	OH	44111	216.671.0577	NDXSalem.com
NDX Stern Empire	1805 W 34th Street	Houston	TX	77018	713.688.1301	NDXSternEmpire.com
NDX Thayer	131 Old Schoolhouse Lane	Mechanicsburg	PA	17055	717.697.6324	NDXThayer.com
NDX Thoele	2801 W St. Germain	St. Cloud	MN	56301	320.252.2070	NDXThoele.com
NDX Twin Cities	1740 Prior Avenue	Falcon Heights	MN	55113	763.786.9121	NDXTwinCities.com

Aesthetic / Specialty

NDX D.H. Baker	2531 Aero Park Drive	Traverse City	MI	49686	800.946.8880	NDXDBaker.com
NDX nSequence	6980 Sierra Center Parkway, Suite 100	Reno	NV	89511	775.827.6645	nSequence.com
NDX Perry & Young	14100 East Evans Avenue	Aurora	CO	80014	303.671.0903	NDXPerryYoung.com

Offshore

Innovative Dental Arts	100 Colony Drive	Irwin	PA	15642	866.305.5434	idarts.net
NDX Pacific	1601 N. Sepulveda Boulevard #628	Manhattan Beach	CA	90266	866.469.0305	NDXPacific.com



NDX EXCLUSIVE PRODUCTS



ClearDream™

ClearDream™

- ClearDream is intended to reduce nighttime snoring and mild to moderate Obstructive Sleep Apnea (OSA) in adults
- Durable: Made with Crystal Clear 450®
- Design: Anterior opening for patient comfort, better air flow and additional tongue space
- Retention: No ball clasps required—the combination of the design and NDX proprietary Clear 450 material achieves sufficient retention



Verotek™
LIFETIME GUARANTEE

Verotek™ Layered

- Porcelain fused-to-zirconia, providing the most aesthetic, natural-looking result
- Hand-layered for internal characterization and suitable for unlimited span

Verotek Aesthetic

- Full-contour zirconia with natural color & translucency
- Ideal for anterior & posterior single units and anterior 3-unit bridges

Verotek FCZ

- Superior strength monolithic zirconia, ideal for bruxers & grinders and suitable for unlimited span



NTI^{TSS}
YOUR MIGRAINE GUARD

The NTI-tss Plus™ can:

- Treat medically diagnosed migraine pain
- Reduce muscle tension
- Prevent bruxism and TMJ syndrome

NDX® NATIONAL DENTEX LABS



RESTORATION GUIDE

Crown & Bridge

	Indications	Characteristics	Contraindications	Material Composition	Preparation Requirements	Cementation
Verotek™ FCZ 	- Posterior Crowns - Bridges (unlimited span) - Suitable for implants	- Metal-free - Superior strength (>1100 mPa) - Good aesthetics	Situations where there is less than 1.0mm of inner occlusal space	- Monolithic medical-grade zirconia - Biocompatible	- Any margin preparation acceptable; slight chamfer preferred 0.8mm axial wall reduction 1.0mm of occlusal reduction	- Remove saliva contamination thoroughly - Use RRG1 or self-adhesive resin cement for maximum retention
Verotek Aesthetic 	- Posterior and anterior crowns - Anterior 3-unit bridges - Suitable for implants	- Metal-free - High strength (780 mPa) - High aesthetics	- Situations where there is less than 1.0mm of inner occlusal space - Bridges or splinted units only up to 3 units	- Monolithic medical-grade zirconia - Generation II Translucency - Biocompatible	- Slight chamfer margin preferred 1.2mm axial wall reduction 1.25mm occlusal reduction preferred	- Remove saliva contamination thoroughly - Use RRG1 or self-adhesive resin cement for maximum retention
Verotek Layered 	- Posterior and anterior crowns - Bridges (unlimited span) - Suitable for implants	- Metal-free - High strength substructure - Premium porcelain overlay	- Situations where there is less than 1.5mm of inner occlusal space - For Verotek Pressed, bridges or splinted units with linear spans > 40mm	- Medical-grade zirconia substructure - Feldspathic porcelain overlay - Biocompatible	- Slight chamfer margin preferred 1.5mm axial wall reduction 1.5mm occlusal reduction	- Remove saliva contamination thoroughly - Use RRG1 or self-adhesive resin cement for maximum retention
PearlPress™ 	- Posterior and anterior crowns - Veneers - Inlays and onlays	- Metal-free - High strength (400 mPa) - Superior aesthetics	- Posterior, 4+ unit bridges or inlay-retained bridges - Very deep subgingival preparations - Patients with substantially reduced residual dentition - Bruxism	Lithium disilicate	1.5mm occlusal and axial wall reduction 1.0mm chamfer margin required	- Bonding: Crowns, veneers, inlays and onlays when tooth preps are primarily supragingival and isolation is possible - Cementation: Indicated for full coverage crowns only
PFM <i>porcelain fused-to-metal</i> 	- Posterior crowns - Posterior bridges	- Moderate strength - Moderate aesthetics	Highly aesthetic cases with minimal prep reduction	Au, Pd, or Ni/Cr alloy substructure with feldspathic porcelain	- Any margin design with 1.0-1.5mm reduction for color 1.5-2.0mm reduction at incisal/occlusal and facial	Conventional cementation of your choice
Full-Cast 	- Posterior crowns - Posterior bridges - Inlays and onlays	- High strength - Very low aesthetics - Kind to opposing dentition		Au, Pd, or Ni/Cr	Feather-edge margin preparations; however any margin preparation may be used	Conventional cementation of your choice

 Posterior Crown
  Inlay/Onlay
  Bridge
  Anterior Crown
  Veneer

RESTORATION GUIDE

Removable Restoration Guide

	Indications	Contraindications	Material Composition	Reline Chairside	Ease of Repair	Overview
Duraflex™	- Flexible/easy insertion - Design versatility - Acrylic allergies - Metal sensitivity	- Need for rigid appliance - Less than 3mm vertical clearance - Flat ridge	Thermoplastic	No	Difficult	The translucent material allows natural coloring of tissue to show through. Material will not gum up on burs while adjusting.
Valplast	- Flexible/easy insertion - Design versatility - Acrylic allergies - Metal sensitivity	- Need for rigid appliance - Less than 5mm vertical clearance - Flat ridge	Nylon Thermoplastic	No	Difficult	Thinner and more translucent. Minimal to no prep required, easy insertion and metal-free comfort.
Framework Partial	- Long and free end saddles - Flat ridge - Few remaining teeth - Strong and stable	- Visible/aesthetics (metal only) - Metal or acrylic allergies (metal only)	Acrylic to Chrome Cobalt or Acrylic to Polymer	Yes	Easier	Can be designed traditionally or with hidden clasps.
Combo	- Flexible desired, but strength needed - Long and free end saddles - Strong and stable - Flat ridge	Lack of vertical clearance	Flexible to Chrome Cobalt	No	Difficult	Utilizes the strength of a cast framework and aesthetics of a flexible. Tooth and tissue supported.
Acrylic Flipper	Temporary treatment partial	- Not long-term - Acrylic allergy	Acrylic	Yes	Easier	Typically used as transitional appliances. Ideal for healing implant placement or extraction.

Denture




Duraflex



Valplast



Framework Partial

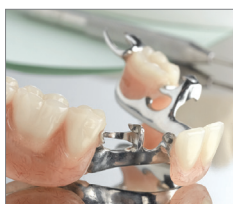

Combo



Acrylic Flipper



RESTORATION GUIDE

Appliances: Therapeutic Migraine Relief

	Indications	Contraindications	What to Send	Overview	Materials
NTI-tss Plus™	- Migraine/tension headaches - Tooth wear from bruxing and clenching - Muscle pain associated with parafunction - Diagnostic treatment planning - TMD	- Advanced periodontal disease - Edge-to-edge bite	- Full arch PVS impressions, models or digital files - Maximum protrusive measurement in mm - Centric occlusion bite	Provides instant and continual incisal guidance in centric clench and any excursions. The standard design only covers the centrals and laterals, but can be extended to cover any number of teeth.	Clear 450® thermoplastic material, providing a perfect fit with no shrinkage and won't absorb stains or odors
NTI-tss Plus Soft	Same as NTI-tss Plus when patient has: - Crowded anteriors - Veneers - Slight mobility	- Advanced periodontal disease - Severely worn centrals - Incompatible occlusal scheme - Severe bruxism	- Full arch PVS impressions, models or digital files - Maximum protrusive measurement in mm - Centric occlusion bite	The NTI-tss Plus Soft features a hard exterior with a soft interior extended to the 2nd bicuspids. It reduces chairtime and increases patient comfort while protecting teeth, muscles and joints from parafunction's devastating effects.	NDX hard/soft thermoformed material with a Clear 450 Discluding Element (DE) affixed

Appliances: Bite Guards

	Indications	Contraindications	What to Send	Overview
Comfort H/S™	- Prevents damage from parafunction - Protect teeth and restorations - Diagnose occlusal concerns	- Short clinical crown - Lack of undercut for retention - Advanced periodontal disease - Multiple edentulous areas	- Full arch PVS impressions, models or digital files 2mm open bite registration	Protect your dentistry and increase patient compliance with the comfort and durability of Comfort H/S. Its hard exterior will withstand your patients' parafunction while the soft liner and economical price make it a prime option.
ThermoFit®	- Protection of natural teeth and restorations - Diagnose occlusal concerns - Prevent damage from parafunction	- Advanced periodontal disease - Mobile teeth	- Full arch PVS impressions, models or digital files 2mm open bite registration	Protect your dentistry and minimize chairtime with this rigid, yet flexible splint. Immerse in hot water for 10 seconds and it drops right in. No internal adjustments needed. When it's in the mouth, it quickly returns to a rigid state.
Crystal Clear®	- Prevents damage from parafunction - Protect teeth and restorations - Diagnose occlusal concerns	Advanced periodontal disease	- Full arch PVS impressions, models or digital files 2mm open bite registration	Made with Clear 450 material. Available in multiple designs and will not absorb stains or odors.

NTI-tss Plus



NTI-tss Plus Soft



Comfort H/S



ThermoFit



Crystal Clear



RESTORATION GUIDE

Appliances: Snoring & Sleep Apnea

	Indications	Contraindications	What to Send	Overview
ClearDream™	ClearDream is intended to reduce nighttime snoring and mild to moderate Obstructive Sleep Apnea (OSA) in adults	- Central Sleep Apnea - Severe respiratory disorders - Loose teeth or advanced periodontal disease - Under 18 years of age - Missing 1st and 2nd molars	- Full arch VPS impressions, models or digital files - A bite registration (The George Gauge bite registration technique is recommended.)	A dorsal appliance that uses 'fins' on the lower arch to advance the mandible. ClearDream maintains an open airway using titratable posterior hardware in the upper arch to adjust the degree the lower jaw is moved forward. Made with Clear 450 material and one piece lower arch construction for better durability.
dreamTAP	- Anti-snoring - Obstructive Sleep Apnea - CPAP intolerance	- Loose teeth or dentures - Central Sleep Apnea	- Full arch VPS impressions, models or digital files - Maximum protrusive measurement in mm - Centric occlusion bite	Made of two custom interocclusal parts that fit over the mandible and maxilla. The lower portion of the device hooks into the upper portion, enabling the mandible to be pulled forward. Titratable using a dial mechanism. Allows greater range of lateral movement for patient comfort.
TAP® 3	- Anti-snoring - Obstructive Sleep Apnea - CPAP intolerance	- Loose teeth or dentures - Central Sleep Apnea	- Full arch VPS impressions, models or digital files - Maximum protrusive measurement in mm - Centric occlusion bite	Made of two custom interocclusal parts that fit over the mandible and maxilla. The upper portion of the device hooks into the lower portion, enabling the mandible to be pulled forward. Titratable using a dial mechanism. Uses a hook mechanism to hold jaw in place.
EMA®	- Snoring - Obstructive Sleep Apnea - CPAP intolerance	- Severe bruxing - Edentulism - Central Sleep Apnea - Mobile teeth	- Full arch VPS impressions, models or digital files - Centric occlusion bite	The primary treatment mechanism of opening the bite and gently moving the mandible forward is achieved using interchangeable elastic straps that offer varying degrees of mandibular advancement. By holding the mandible forward, the EMA maintains a clear airway to reduce snoring and improve breathing.

ClearDream



dreamTAP



TAP 3



EMA



EDUCATIONAL OPPORTUNITIES



CONTINUING EDUCATION PROGRAMS

National Dentex provides Continuing Education opportunities to thousands of clinicians around the country on an annual basis. We offer a variety of seminars and webinars on relevant topics in dentistry.

Lunch & Learn Opportunities

We offer many CE programs for our clients and their staff through our Lunch & Learn programs conveniently presented in your office. We'll provide lunch and you earn CE credit.

Video Library

At National Dentex, we are committed to providing our clients and their patients access to useful and relevant content in regards to dental products and services that we offer.

We invite you to explore our ever-growing list of online videos to help inform you of trends in the industry and opportunities to grow your practice.

Document Library

NDX offers reference materials for your staff and office. Review our content and check back for new titles.

- Helpful 'How-To' Guides
- Troubleshooting Guides
- Product Guides
- Informative Product-Related Articles

Contact us today at ndxmarketing@nationaldentex.com or visit nationaldentex.com/education to learn how we can provide training for you and your staff.

DIGITAL DENTISTRY

iTero® ■ TRIOS® ■ Carestream ■ Cerec ■ 3M™ ■ Planmeca®



6 SYSTEMS. 1 LABORATORY NETWORK. DIGITAL IMPRESSIONS MADE EASY.

Our NDX consultants can help guide you through the digital evolution with dental programs, data plans and equipment for your practice. National Dentex Laboratories not only process digital data from all major scanning systems, making the transition to digital dentistry easy for your practice, but we also provide consulting services to help you choose the scanning system and programs that work best for your practice. Some of the benefits of moving to digital dentistry include:

Remake Reduction

Digital impressions provide increased accuracy and the ability to identify issues ahead of time, which reduce potential fit and margin discrepancies.

Laboratory Relationship

Enjoy working with a lab with the ability to handle your toughest cases and any questions regarding digital dentistry.

Flexibility

Choose the digital impression system that works best for your office. National Dentex Laboratories accept scans from all major brands including 3M True Definition, Cerec, iTero, 3Shape TRIOS, Planmeca/ E4D Technologies and Carestream.

Expedited Turnaround

Digital dentistry delivers your case to the lab faster than traditional methods, which reduces the turnaround time necessary for each case.

Cost Effective

Save money on materials typically used during the traditional impression-taking process. Most model-less cases allow for reduced lab fees.

Increased Patient Comfort

Digital impressions are much quicker and therefore more comfortable than the traditional impression-taking process as there is no need for excess impression material.

INTRAORAL SCANNER INFO CHART

		Cerec (Omnicam or Bluecam)	iTero®	True Definition (COS/3M)	Planmeca Planscan® (formerly E4D)	TRIOS® (3Shape)	Carestream
Process	Add Your NDX Lab	Create a Cerec Connect account, then place your NDX Lab in "My Favorite labs" folder	Add account ID# to scanner menu or contact Cadent/ iTero tech support	Create a 3M Connection Center account and call 3M Digital Service to add your NDX Lab to the drop-down menu	In the CAD/CAM module of the active patient, select the desired case file, then click DDX Export	Contact lab for its email to add it to your 3Shape Communicate Account	For the first case, contact your local NDX lab for the appropriate e-mail address
	Send a File	Select your NDX Lab from your favorites folder	Select your NDX Lab from your drop-down menu	Select your NDX Lab from your list of labs	Send the files through the Henry Schein DDX platform	Select your NDX Lab from your list of available labs	Select your NDX Lab from the drop-down menu in CS Connect
	Request Model-less	Specify the case is model-less in the case notes section		Ask True Def to add your NDX Lab model-less option to your drop-down lab menu	Specify the case is model-less in the case notes section		Specify that the case is model-less in the comments box
Available Products	Via Scan Submission	All Ceramics, PFM, Implants, Verotek, Lithium Disilicate Monolithic Restorations and Select Removable Products	All Ceramics, PFM, Implants, Nightguards, Verotek, Lithium Disilicate Monolithic Restorations and Select Removable Products		All Ceramics, PFM, Implants, Verotek, Lithium Disilicate Monolithic Restorations and Select Removable Products		
	Require Models	Implants, Nightguards and Removables					
	Can be made model-less	Verotek and Lithium Disilicate Monolithic Restorations					
Warranty	Models	Standard Warranty					
	Model-less	Standard Warranty					
Do you keep my scan on file?		60 Days					

Contact your local NDX lab to find out about NDX model-less programs.

EXPECT MORE.

CHECK OUT ALL WE HAVE TO OFFER

SERVICES

- ✓ Personalized Customer Service
- ✓ On-Site Support
- ✓ Mail Order Delivery
- ✓ Continuing Education
- ✓ Resource Libraries
- ✓ Offshore Solutions
- ✓ Practice Consulting
- ✓ Business Seminars

PRODUCTS

- ✓ Crowns and Bridges
- ✓ Dentures and Partials
- ✓ Specialty Appliances
- ✓ Migraine Therapy
- ✓ Implant Restorations
- ✓ Dental Sleep Medicine
- ✓ Digital Surgical Guides
- ✓ Orthodontics



SERVICES

Expect NDX services to revolutionize your practice. We deliver consistent quality and up-to-date expertise that support your goal of delivering superior care and will broaden the scope of your business. Learn more about our unique offerings from our expansive list of services.

Dedicated Customer Service

NDX provides timely pick-up and delivery and status of cases

Chairside Expertise and Shade Taking

Expert custom shading and denture conversions

On-Site Certified Technicians and Technical Support Managers

Available to assist you on your most challenging cases

Mail Order Delivery

NDX provides quality restorations delivered anywhere in the United States

State-of-the-Art Milling Centers

Design your cases for accurate restorations

Comprehensive Digital Workflow

Our labs fabricate your restorations and products utilizing the most up-to-date CAD/CAM, milling and 3D printing technologies

Intraoral Scanners

NDX labs accept IOS scans from all major scanners on the market, saving your office time and money

Restoration Guide Library

NDX can furnish easy-to-use guides for product selection

Dental Education Video Library

Offers educational videos on procedures, products, business practices and more

Continuing Education Programs

On-demand webinars and local seminars to keep you up-to-date and your practice growing

Lunch and Learns

Continuing education opportunities in the comfort of your office

Hands-On Workshops

Get hands-on experience with the latest technology and learn best practices

Territory Managers

Serve as your personal lab resource

Emodel Orthodontic Services

Save space in your office by storing models digitally

Domestic and Offshore Solutions

Offer options to fit your quality and cost requirements

Visit NationalDentex.com/EXPECTMORE or reach out to your local NDX lab or Territory Manager for additional information. *Not all programs available at all locations.*



PRODUCTS

NDX provides products that exceed the expectations of your patients. Grow your practice with our exclusive offerings such as nSequence Guided Prosthetics® and the NTL-tss Plus® Migraine Guard. Find the latest options for flexible partials and a variety of custom appliances to treat sleep apnea. We offer products with a lifetime guarantee such as Verotek™ and have strong relationships with implant companies that span many years. With NDX products, expectations meet reality.

Dentures and Partial

- Flexible Partial
- Cast Partial
- ClearFrame™
- Full Dentures
- Ultraire AKP™

Specialty Appliances

- Comfort H/S™ Bite Splint
- Crystal Clear®
- Thermo-Guard™

Dental Sleep Medicine

- ClearDream®
- dreamTAP®
- TAP® 3
- EMA®

Crowns and Bridges

- Verotek Full Contour Zirconia*
- Verotek Layered Zirconia*
- Verotek Aesthetic Zirconia*
- PearlPress™
- IPS e.max®
- Porcelain-Fused-to-Metal
- Full-Cast

**All Verotek options come with a lifetime guarantee.*

Orthodontics

NDX offers complete orthodontic capabilities

Migraine Therapy

The NTL-tss Plus is the only FDA-cleared dental appliance for the treatment of migraine pain

Implants

NDX has experience working with all implant systems from a variety of manufacturers and can offer OEM for cases where original manufacturer parts are requested for a single or a multi-unit implant case

Digital Surgical Guides

nSequence Guided Prosthetics, our advanced solution for implant-supported full-arch restorations, is a comprehensive guided stackable system that saves substantial chair time



EXPECT...

- Simplicity from working with one partner who offers all the products and services you need
- Personalized service and support from the NDX national network of labs
- Superior aesthetic products and the most advanced technology in digital dentistry and surgical planning
- Educational opportunities that continually enhance your professional and business knowledge
- Growth from a range of services that increase patient volume, practice efficiencies and profitability

EXPECT MORE when partnering with NDX.

To learn more, please call 844.295.6305 or visit NationalDentex.com/EXPECTMORE



11601 Kew Gardens Avenue, Suite 200

Palm Beach Gardens, FL 33410

NationalDentex.com

To learn more, please call 844.295.6305
or visit NationalDentex.com/EXPECTMORE

The Power of Our Network

National Dentex Labs is the largest network of dental laboratories in the nation. Together, we offer a range of products, services and technologies to address any dental restoration. To explore our full range of products, please visit: www.NationalDentex.com.



We are dedicated to building personal relationships with our services: local representation, drivers, chairside expertise and continuing education.



Offering simplicity, value and full capabilities, our labs ship anywhere in the USA—a perfect solution for both DSO's and individual practices.



When aesthetics matter most, our most experienced craftsmen provide meticulously prepared, beautiful dental restorations.



NDX offers cutting-edge surgical guides and digital preparation services, such as nSequence® Guided Prosthetics®.

KENTUCKY - Louisville
NDX[®]  Affinity

PENNSYLVANIA - Irwin
NDX[®]  Albenzi

SOUTH DAKOTA - Mitchell
NDX[®]  Bauer

ARIZONA - Phoenix
NDX[®]  Continental

MICHIGAN - Lansing
NDX[®]  Dental Art

MINNESOTA - Farmington
NDX[®]  Dexterity

MICHIGAN - Traverse City
NDX[®]  D.H.Baker

DELAWARE - New Castle
NDX[®]  Dodd

FLORIDA - Orlando
NDX[®]  Florida

ARKANSAS - Heber Springs
NDX[®]  Green

NEW HAMPSHIRE - Manchester
NDX[®]  H&O

NEW MEXICO - Albuquerque
NDX[®]  Ideal

MISSOURI - Fenton
NDX[®]  Keller

NEBRASKA - Omaha
NDX[®]  Kiess Kraft

WISCONSIN - De Pere
NDX[®]  Lords

INDIANA - Indianapolis
NDX[®]  Lumident

MISSOURI - Lee's Summit
NDX[®]  Mallow-Tru

NEVADA - Reno
NDX[®]  nSequence

IOWA - Dubuque
NDX[®]  Oral Arts Iowa

CALIFORNIA - Manhattan Beach
NDX[®]  Pacific

COLORADO - Aurora
NDX[®]  Perry & Young

TENNESSEE - Smyrna
NDX[®]  Peterman

TENNESSEE - Athens
NDX[®]  Rogers

WISCONSIN - Waukesha
NDX[®]  Saber

OHIO - Cleveland
NDX[®]  Salem

TEXAS - Houston
NDX[®]  Stern Empire

PENNSYLVANIA - Mechanicsburg
NDX[®]  Thayer

MINNESOTA - St. Cloud
NDX[®]  Thoele

MINNESOTA - Falcon Heights
NDX[®]  Twin Cities



The Network of Artisans and Innovators

National Dentex Labs
11601 Kew Gardens, Suite 200
Palm Beach Gardens, FL 33410

www.NationalDentex.com